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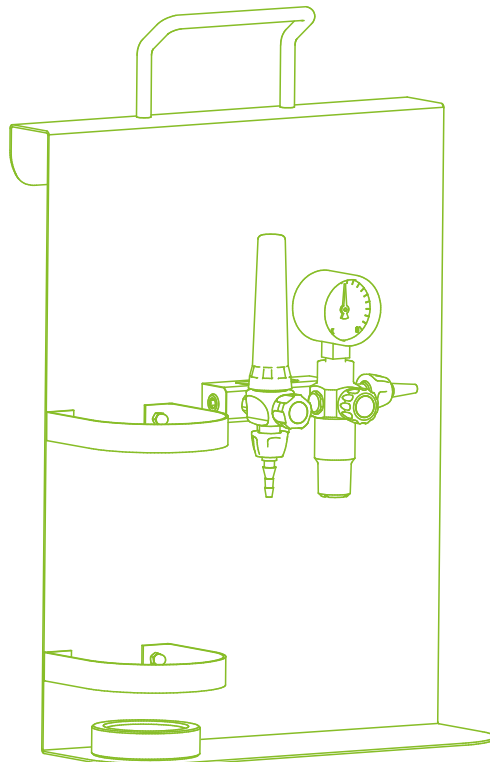
Habicht emergency supply unit

Dear customers,

thank you for purchasing this Greggersen product. If you have any questions or require information, please contact our sales and support team.

+49-(0)40 739 357-0, sales@greggersen.de

Preface

Please read the instructions for use thoroughly before using the medical device. Only use it once you have completely understood the application and function!

Intended use

For insufflation and extraction if central oxygen or medical air supply is not possible (emergencies, patient transfer, home-care, etc.).

Necessary accessories

- variant with Kolibri 2 or 3 l O₂ cylinder
- variant with tube flow meter 3 l O₂ cylinder
- pressure reducer
- low pressure hose

Operation

Connect the pressure reducer to the O₂ bottle and attach it to the emergency supply unit at the NIST connection using a low-pressure hose line. The flow meter and the suction unit can now be put into operation. Always close the cylinder valve during longer interruptions. Remove pressure from the unit when replacing the cylinder, detach the hose from the pressure reducer and unscrew the connection nut of the pressure reducer from the cylinder valve by hand (no tools).

Cleaning

It is recommended to clean the device by wiping disinfection after each use in accordance with the cleaning agent manufacturer's instructions. Commercially available disinfectants that are compatible with polycarbonate should be used. When cleaning, make sure that no liquid gets into the appliance. The device must be fully dried prior to reuse. Do not soak the device with disinfectant! The device may be damaged. When using the Greggersen hose connector, it must be cleaned and disinfected in the same manner as the device. The hygiene plan of the institution must be observed. The device is designed for multiple use with different patients.

Warnings

- Prior to the use of the device, the user must be familiar with the product and the instructions for use and must have understood them.
- The instructions for use of the flow meters, suction units and pressure reducers apply.
- The emergency supply unit Habicht with tube flow meter is not suitable for use with 2 liter bottles.
- Simultaneous use of both functions influences the performance and measurement accuracy.
- Never use tools when connecting the unit and the accessories; all connections can be tightened or loosened by hand. Always unscrew accessories when the pressure reducer is pressureless.
- The flowmeter will only display correctly if the measuring tube is vertical.
- Due to the oxidising effect of oxygen, special care must be taken when using oxygen-powered suction units.

Technical data

Type: Portable stainless steel frame, for suspension from an appliance rail or patient bed; with holder for gas cylinder of max. 3 liter (max. Ø = 100mm); Consisting of ejector regulator and tube flowmeter or flowmeter Kolibri


Dimensions: approx. 500x330x150mm (HxWxD)

Medical Products Directive 93/42/EEC; Classification: IIa

Maintenance and inspection

The emergency supply unit must be serviced at least every 5 years. An inspection with a functional test and a leak test has to be performed at least every 2 years. In case of unusually heavy-duty use, shorter servicing and inspection intervals may be required.

Servicing and repairs must only be carried out by authorised expert personnel. Only use genuine Greggersen spare parts.

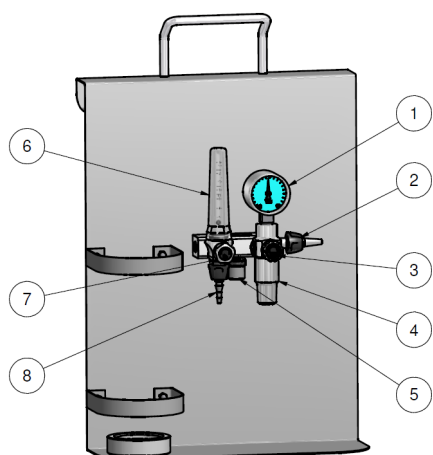1. Gauge
2. Vacuum inlet (9/16" UNF)
3. Needle valve vacuum
4. Ejector regulator
5. Pressure inlet (NIST)
6. Flowmeter tube / Kolibri
7. Needle valve flowmeter
8. Flow outlet (9/16" UNF)