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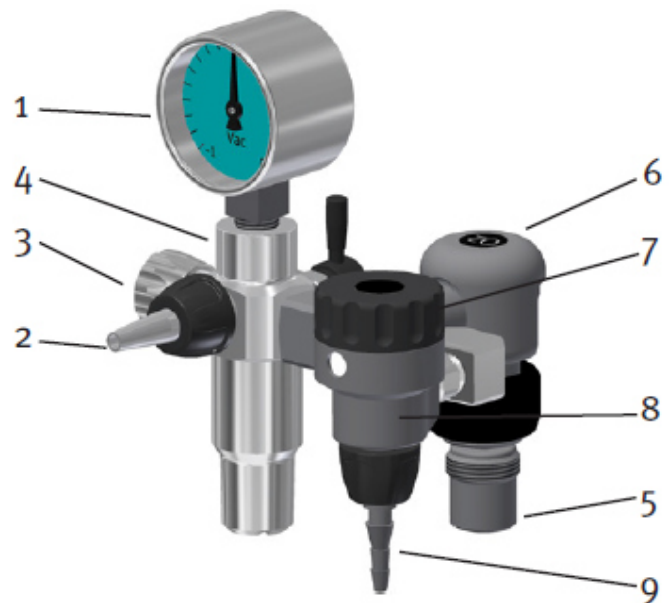
Habicht rail unit

Dear customers,

thank you for purchasing this Greggersen product. If you have any questions or require information, please contact our sales and support team.

+49-(0)40 739 357-0, sales@greggersen.de

Preface

Please read the instructions for use thoroughly before using the medical device. Only use it once you have completely understood the application and function!

Intended use

For insufflation and extraction.

Necessary accessories

- humidifiers
- atomizers
- one-way humidifiers
- overflow bottle
- flowmeter nozzle
- vacuum nozzle
- low-pressure hose system

Operation

Fasten the rail unit to the equipment rail. Connect the device to a gas outlet using a NIST hose assembly. Close the metering valve on the flowmeter and the ejector regulator. Connect the humidifier bottle, overflow bottle and any other accessories.

Use the handwheel to set the desired gas flow. During use, ensure that there is a sufficient water level in the humidifier or nebulizer. Set the ejector regulator to the desired suction level. For short interruptions, close the valves. For longer interruptions, always disconnect the hose assembly from the gas outlet.

Cleaning

It is recommended to clean the device by wiping disinfection after each use in accordance with the cleaning agent manufacturer's instructions. Commercially available disinfectants should be used. When cleaning, make sure that no liquid gets into the appliance. The device must be fully dried prior to reuse. Do not soak the device with disinfectant! The device may be damaged. When using the Greggersen hose connector, it must be cleaned and disinfected in the same manner as the device. The hygiene plan of the institution must be observed. The device is designed for multiple use with different patients.

Warnings

- The device and all connections must be free from oil, grease and lubricants, as otherwise there may be a risk of fire or explosion.
- Simultaneous use of several functions may affect performance and measurement accuracy. Simultaneous use of the flowmeter, suction unit and gas outlet may result in performance below the specified values.
- Never use tools when connecting accessories. All connections can be tightened or loosened by hand. Only connect or disconnect accessories when the system is depressurized.
- Do not forcibly close the fine metering valves, as this may damage the valve seat.
- Subject to technical changes; other versions are available on request.
- The owner/operator of the device must ensure that all users of the device are familiar with the instructions for use and the operation of the device.
- Warranty claims are excluded in the event of damage or malfunction caused by improper transport, improper storage or use not in accordance with the intended purpose.
- Fluctuations in inlet or back pressure, as well as ambient conditions, may affect the flow rate.
- Observe the instructions for use of the flowmeter and ejector regulator.
- Due to the combustion-supporting properties of oxygen, special care must be taken when using oxygen-powered suction units.

Technical data

Type:	Rail unit with flowmeter and ejector regulator; all units of chrome-plated, all-metal construction; consisting of: - flowmeter kolibri, latch-down flow rate adjustment for medical gases, indication via flow figures - ejector regulator with fine-control valve and vacuum gauge
Inlet:	NIST according to DIN EN ISO 5359
Capacity:	Flowmeter adjustable flow rate in steps between 0 ... 15 l/min Ejector regulator approx. -0.85 bar, suction power max. 15 l/min., oxygen consumption max. 30 l/min.
Guidelines:	flowmeter DIN EN ISO 15002 ejector DIN EN ISO 10079-3 NIST-connection DIN EN ISO 5359
Medical Products Directive 93/42/EEC; Classification: IIa	

Maintenance and inspection

External disinfection with a disposable cloth is sufficient for all chrome-plated metal parts. Clean plastic parts with a common disinfectant solution and follow manufacturers instructions. Let equipment dry completely before using it again.

Prior to each use, make a visual check (external damage). A function and leak test must be performed at least once a year. Maintenance must be done every five years (replacement of wear parts). Failures and damages may only be remedied by authorized expert personnel. Do repairs only on the pressureless unit. Use only original **GREGGERSEN** spare parts.

1. vacuum gauge
2. vacuum nozzle
3. ejector regulator
4. needle valve
5. NIST-connection
6. rail clamp
7. handwheel
8. flowmeter Kolibri
9. flowmeter nozzle