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
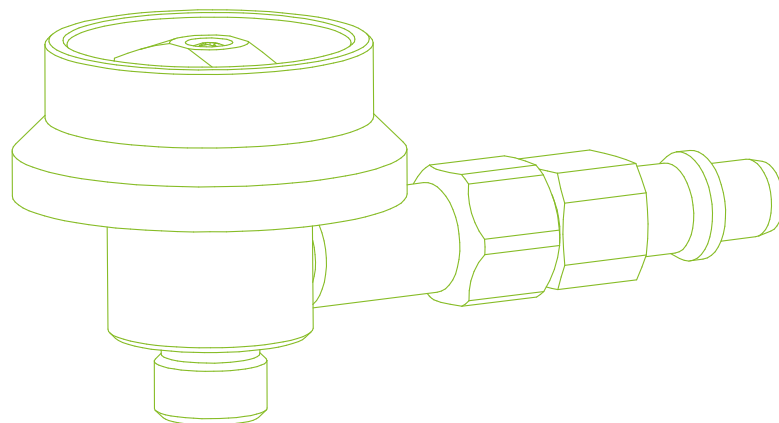
Medical atomizer valve

Dear customers,

thank you for purchasing this Greggersen product. If you have any questions or require information, please contact our sales and support team.

+49-(0)40 739 357-0, sales@greggersen.de

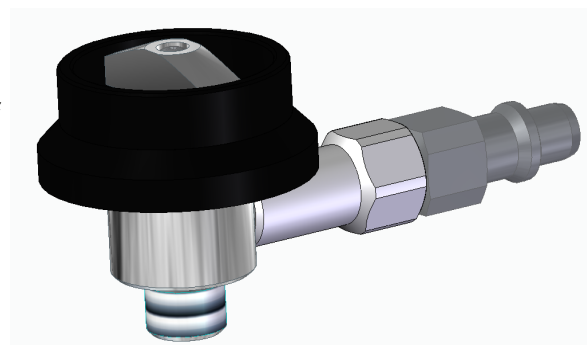Preface

The Medical atomizer valves from Greggersen Gasetechnik GmbH comply with the regulations of the Medical Device Directive 93/42/EEC and are classified as active, non-invasive devices of Class IIa.

Operation

Close the quick shut-off valve and connect the plug-in valve directly to the gas outlet with the outlet connection pointing downwards. Connect the Inhalette to the outlet connection and fill it with the liquid aerosol solution. For inhalation, hold the Inhalette in your hand, fully open the quick shut-off valve and begin inhalation.

After use, close the quick shut-off valve and place the Inhalette back onto the holder of the plug-in valve. The plug-in valve may remain connected to the gas outlet, preferably in the parking position.

When the Inhalette is not in use, it can be placed on the holder of the plug-in valve.

Warnings

- Always keep all connections free of oil and grease!
- Do not use force to close the quick shut-off valve; the valve seat may be damaged.
- The metering screw is matched to the inlet pressure, type of gas and required flow rate. Using other metering screws will affect the accuracy of the flow rate!
- Never use tools when connecting accessories. All connections can be tightened or loosened by hand. Accessories must only be connected or disconnected when the system is depressurized.
- The medical atomizer valve may only be used by trained personnel.
- Subject to technical modifications; other versions are available on request.
- Fluctuations in inlet or back pressure as well as changes in ambient conditions may affect the flow rate.
- Warranty claims are excluded in the event of damage or malfunctions caused by improper transport, improper storage or use other than as intended.

Cleaning

It is recommended to clean the device by wiping disinfection after each use in accordance with the cleaning agent manufacturer's instructions. Commercially available disinfectants should be used. When cleaning, make sure that no liquid gets into the appliance. The device must be fully dried prior to reuse. Do not soak the device with disinfectant! The device may be damaged. When using the Greggersen hose connector, it must be cleaned and disinfected in the same manner as the device. The hygiene plan of the institution must be observed. The device is designed for multiple use with different patients.

Maintenance/Service

External wipe disinfection is sufficient for all chrome-plated brass parts. Clean plastic parts with a disinfectant solution. Allow the device to dry completely before using it again.

The medical atomizer valve must be visually inspected for external damage before each use. The devices should be serviced at least every five years (functional and leak test). Shorter maintenance intervals may be required in the event of unusually heavy use. Faults and damage may only be repaired by authorized qualified personnel. Use only original GREGGERSEN spare parts!

Technical data

Design: Chrome-plated all-metal construction
Plug-in device for compressed air
Quick shut-off valve

Flow rate: 5 l/min $\pm$ 20% at an inlet pressure of 450 kPa $\pm$ 50 kPa

Connection: Inlet: Plug-in connectors according to DIN 13260 Part 2
Outlet: 9/16-18 UNF threaded connection

Standards: DIN EN ISO 15002, DIN 13260-2
Medical Device Directive: 93/42/EEC; Classification: IIa

Accessories

900619 Hose connection, 9/16" hose barb – compressed gas