


Automatic Manifold MediControl MC2100E

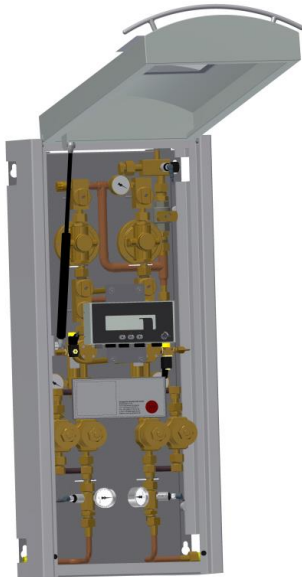

Manufacturer: Greggersen

Product 326100

Number:

Description

Automatic Manifold MediControl 2100E

Flow Rate: 100 m³/h

Inlet Pressure: max. 20000 kPa (200 bar)

Outlet Pressure: 500 kPa +/- 50 kPa (5 bar +/- 0.5 bar)

Electr. Connection: 230 V AC

System Language: German

Mode: 2 Cylinder Banks

Gas Type: ++++++

Consisting of:

1 x MC2100

1 x Pipe Set MC2100

Microprocessor-controlled gas supply system for 2-sided cylinder supply, design in accordance with DIN EN ISO 7396-1, LCD display with operating status indication and service messages in plain text, forwarding of

messages via volt-free contacts, sensor-based pressure monitoring, 2-stage pressure reduction, duplicated 2nd stage, pneum. priority switching in case of power failure, with pressure gauges and emergency supply option (NIST), fully assembled on mounting plate incl. cover