


Automatic Manifold MediControl MC2025R

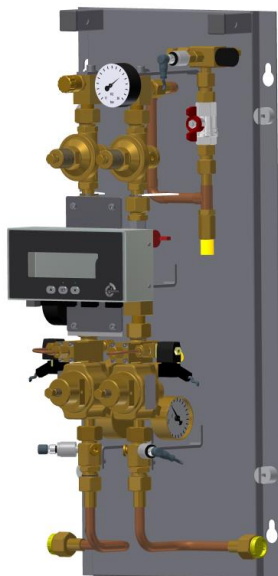

Manufacturer: Greggersen

Product 326026

Number:

Description

Automatic Manifold MediControl 2025R

Flow Rate: 25 m³/h

Inlet pressure: max. 200 bar

Outlet pressure: 5 bar +/- 0.5 bar

Electr. connection: 230 V AC

System language: German

Mode: 2 cylinder banks and reserve supply

Gas type: ++++++

Consisting of:

1 x MC2025

1 x Set of pipes MC2025 + inlet pipe reserve

Microprocessor-controlled gas supply system for 2-sided cylinder supply, design according to DIN EN ISO 7396-1, LCD display with indication of operating status and service messages in plain text, forwarding of

messages via potential-free contacts, sensory pressure monitoring, 2-stage pressure reduction, 2nd stage in dual design, pneum. priority switching in case of power failure, with pressure gauges and emergency feed facility (NIST), fully mounted on mounting plate incl. cover; incl. connection facility for reserve supply