


Vacuum Reg. Pirol 2.0 -90,VAC,RAIL,CH

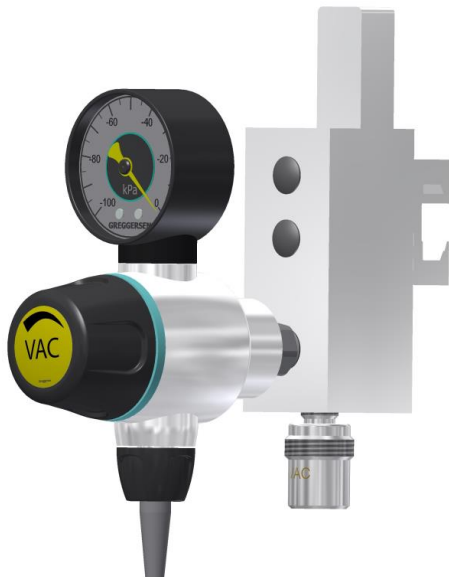

Manufacturer: Greggersen

Product 505641

Number:

Description

Vacuum Regulator Pirol -90, Rail-Mounted Unit CH, NIST
with spring-loaded rail clamp for up to 50 mm rails

Medical vacuum regulator for dosing vacuum.

With 360° rotatable vacuum gauge; large, ergonomic
adjustment handwheel; rail-mounted unit incl. spring-loaded
rail clamp.

The gas inlet is integrated into the rail clamp (NIST).

The adjustment mechanism enables stepless, precise and
intuitive regulation of the vacuum. The rotatable
vacuum gauge allows the user to read the
set value from an optimal viewing angle.

The regulator unit meets the requirements of
DIN EN ISO 10079-3 (high vacuum / high flow).

Material:

Housing: aluminium, anodised

Handwheel: plastic

Rail clamp: aluminium, anodised

Technical data:

Gas type: Vacuum / VAC

Inlet: NIST housing according to DIN EN ISO 18082

Control range: 0 to -90 kPa

Flow rate: >20 l freeflow/min

Outlet: 9/16-18 UNF with vacuum tubing nipple