


## Automatic Manifold MediControl MC2025P

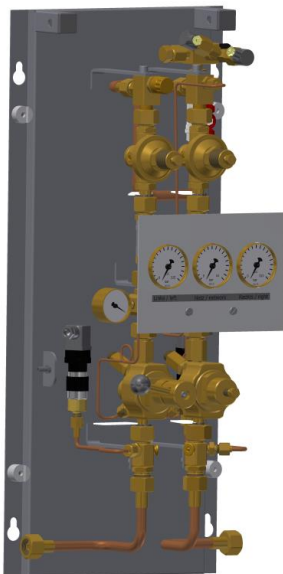

**Manufacturer:** Greggersen

**Product** 327025

**Number:**

### Description

Automatic Manifold MediControl 2025P

Flow rate: 25 m<sup>3</sup>/h

Inlet pressure: max. 200 bar

Outlet pressure: 5 bar +/- 0.5 bar

Gas type: ++++++

Consisting of:

1 x MC2025P

1 x Pipe set MC2025P

Pneumatically controlled gas supply system for two-sided cylinder bank supply, design according to DIN EN ISO 7396-1, automatic changeover, priority control with lever, pressure gauge indication for both supply sources and line pressure, pressure monitoring with contact gauges, two-stage pressure reduction, second stage duplicated, with control pressure gauges, high-pressure and low-pressure contact gauges and emergency feed facility (NIST,) completely mounted on mounting plate, housing with Viewing Window