


Automatic Manifold MediControl MC2100 without Tank

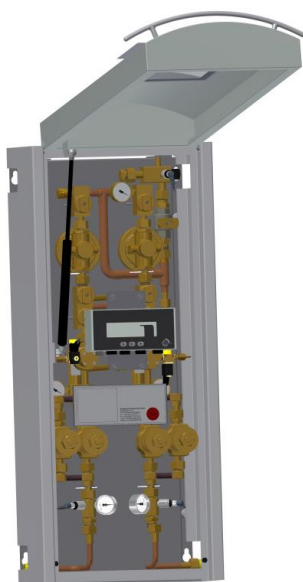

Manufacturer: Greggersen

Product 903700

Number:

Description

Automatic Manifold MediControl 2100T

Flow rate: 100 m³/h

Inlet pressure: max. 20,000 kPa (200 bar)

Outlet pressure: 500 kPa +/- 50 kPa (5 bar +/- 0.5 bar)

Elec. connection: 230 V AC

System language: German

Mode: Tank and 2 cylinder banks

Gas type: ++++++

Consisting of:

1 x MC2100

1 x Pipe set MC2100

Microprocessor-controlled gas supply system for 2-sided cylinder supply, design in accordance with DIN EN ISO 7396-1, LCD display with indication of operating status and service messages in plain text, forwarding of messages via potential-free contacts, sensor-based pressure monitoring, 2-stage pressure reduction, 2nd stage redundant, pneum. priority switching in the event of power failure, with pressure gauges and emergency feed option (NIST), fully assembled on mounting plate incl. cover