


## Vacuum Reg. Pirol 2.0 -30,VAC,RAIL,CH

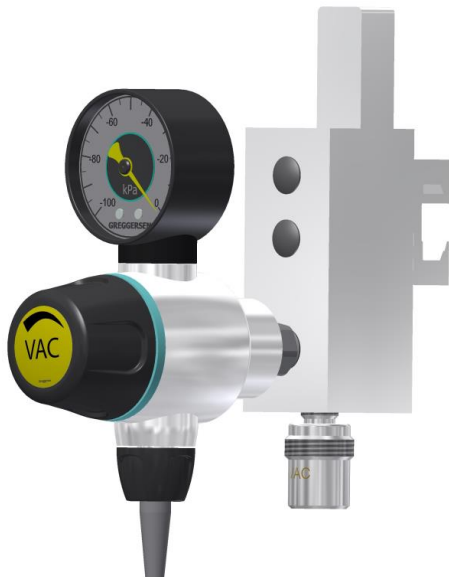

**Manufacturer:** Greggersen

**Product** 505656

**Number:**

### Description

Vacuum Regulator Pirol -30, Rail-Mounted Unit CH, NIST  
with spring-loaded rail clamp for up to 50 mm rails

Medical vacuum regulator for dosing vacuum.

With 360° rotatable vacuum gauge; large, ergonomic  
adjustment handwheel; rail-mounted unit incl. spring-loaded  
rail clamp.

The gas inlet is integrated into the rail clamp (NIST).

The adjustment mechanism enables stepless, precise and  
intuitive vacuum regulation. The rotatable  
vacuum gauge enables the user to read the  
set value from an optimum viewing angle.

The regulator unit meets the requirements of  
DIN EN ISO 10079-3 (high vacuum / high flow).

Material:

Housing: Aluminium, anodised

Handwheel: Plastic

Rail clamp: Aluminium, anodised

Technical data:

Gas type: Vacuum / VAC

Inlet: NIST housing according to DIN EN ISO 18082

Control range: 0 to -30 kPa

Flow rate: >20 l free flow/min

Outlet: 9/16-18 UNF with vacuum tubing nipple