


Automatic Manifold MediControl MC2100R

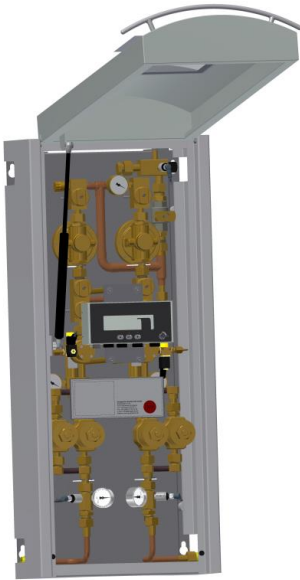

Manufacturer: Greggersen

Product 326102

Number:

Description

Automatic Manifold MediControl 2100R

Flow rate: 100 m³/h

Inlet pressure: max. 20000 kPa (200 bar)

Outlet pressure: 500 kPa +/- 50 kPa (5 bar +/- 0.5 bar)

Electr. connection: 230 V AC

System language: German

Mode: 2 cylinder banks and reserve

Gas type: ++++++

Consisting of:

1 x MC2100

1 x pipe set MC2100 + inlet pipe reserve

Microprocessor-controlled gas supply system for 2-sided cylinder supply, design according to DIN EN ISO 7396-1, LCD display with indication of operating status and service messages in plain text, forwarding of

messages via potential-free contacts, sensor-based pressure monitoring, 2-stage pressure reduction, dual 2nd stage design, pneum. priority switching in the event of power failure, with pressure gauges and emergency supply option (NIST), fully assembled on mounting plate incl. cover, incl. connection facility for reserve supply