


Filter Panel G1 Vacuum

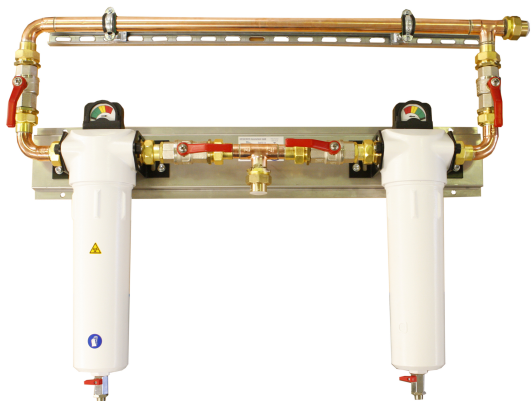

Manufacturer: Greggersen

Product 903510

Number:

Description

Filter Panel Size 1 Vacuum

Duplex bacterial filter stage

up to 65 m³/h vacuum

To protect the vacuum vessel and pumps against contamination by micro-organisms and foreign particles as well as for filtering the drawn-in exhaust air.

For use in a medical vacuum supply system in accordance with DIN EN ISO 7396-1.

Redundant design of the filter stage. This allows filter replacement during maintenance work without interrupting the vacuum supply.

Scope of delivery consists of:

4 ball valves 1"

2 bacterial filters with attached secretion jar

1 wall-mounting plate 1100 x 150mm (WxH) and pipe clamps
for pipe support

Connecting and linking pipes copper pipe 28mm

Technical data:

Capacity: 65 m³/h

Dimensions: 1150 x 900 x 150mm (WxHxD) incl. filter

Connection: Copper pipe Ø 28x1.5mm

Efficiency: 99.99% at 0.01μ

Weight approx.: 20 kg