


Automatic Manifold MediControl MC2025E

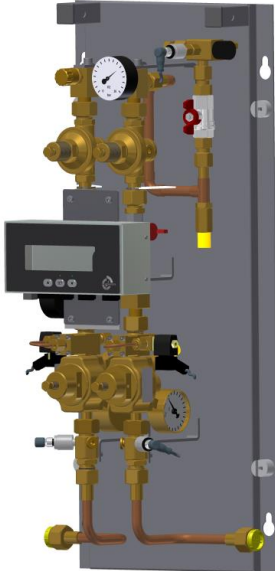

Manufacturer: Greggersen

Product 326025

Number:

Description

Automatic Manifold MediControl 2025E

Flow rate: 25 m³/h

Inlet pressure: max. 200 bar

Outlet pressure: 5 bar +/- 0.5 bar

Elec. connection: 230 V AC

System language: German

Mode: 2 cylinder banks

Gas type: ++++++

Consisting of:

1 x MC2025

1 x Pipe Set MC2025

Microprocessor-controlled gas supply system for 2-sided cylinder supply, design according to DIN EN ISO 7396-1, LCD display with plain text operating status and service messages, forwarding of messages via

potential-free contacts, sensor-based pressure monitoring, 2-stage pressure reduction, 2nd stage duplicated, pneum. priority circuit in case of power failure, with pressure gauges and emergency supply facility (NIST), fully assembled on mounting plate incl. cover