


Ventus 22 mm Valve Insert +NIST AIR7

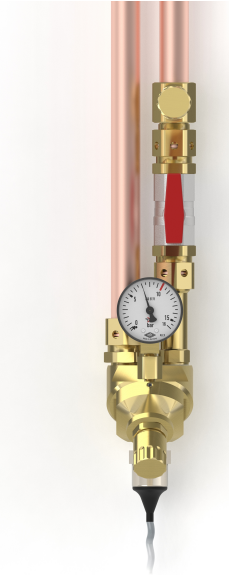

Manufacturer: Greggersen

Product 900798

Number:

Description

Valve insert for Ventus 22 mm

Gas type : AIR7

Design complies with DIN EN ISO 7396-1 and HTM 02-01.

Valve assembly consisting of:

- Shut-off valve DN 20, G3/4"
- Copper pipe 22mm inlet / outlet
- Pressure gauge 0...16 bar
- Pressure sensor 0...16
- NIST emergency inlet point AIR-800 upstream of the shut-off valve
- NIST emergency inlet point AIR-800 downstream of the shut-off valve
- Physical separation point
- Gas type labelling AIR7

Properties

Fulfilled Standard	DIN EN ISO 7396-1 & HTM 02-01
Gas Type	AIR 7
NIST Emergency Inlet Point	2 NIST Emergency Supply Points