


Ventus 22mm 2-way Front Part+Base Part

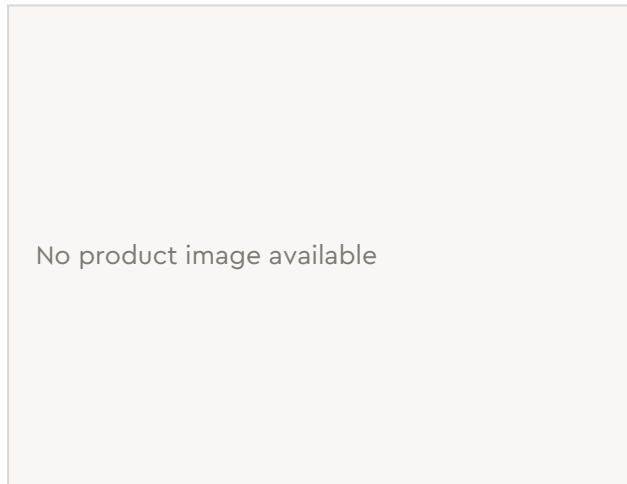

Manufacturer: Greggersen

Product 900878

Number:

Description

Area Control Unit VENTUS 22mm

2-way, Front Part + Base Part

Gas type: *** Please specify ***

Language of Instructions for Use: German, English, Dutch, Danish, Hungarian (please delete inapplicable options!)

Design conforms to DIN EN ISO 7396-1 and HTM 02-01.

Front Part assembly consisting of:

- Sheet steel door, powder-coated
- Side-hinged opening
- 2-way alarm unit
- Magnetic foil to cover the warning
- Plaster compensation up to 25mm
- Power supply unit 230 V AC

Base Part assembly consisting of:

- Installation housing for flush mounting

The alarm unit monitors:

- Pipeline pressure too high
- Pipeline pressure too low

Dimensions

Inlet & Outlet: 22 mm CU pipe

Front Part (WxHxD): 390 x 530 x 9 mm

Base Part (WxHxD): 330 x 470 x 77 mm

Alarm limits of the warning:

O₂, AIR, N₂O, CO₂: 500 kPa ± 20%

VAC: -40 kPa

AIR-800: 800 kPa ± 20%